


Efficient, Versatile, Ergonomic

Epredia Histostar Embedding Center



Epredia HistoStar Embedding Workstation

Cool, comfortable and in control

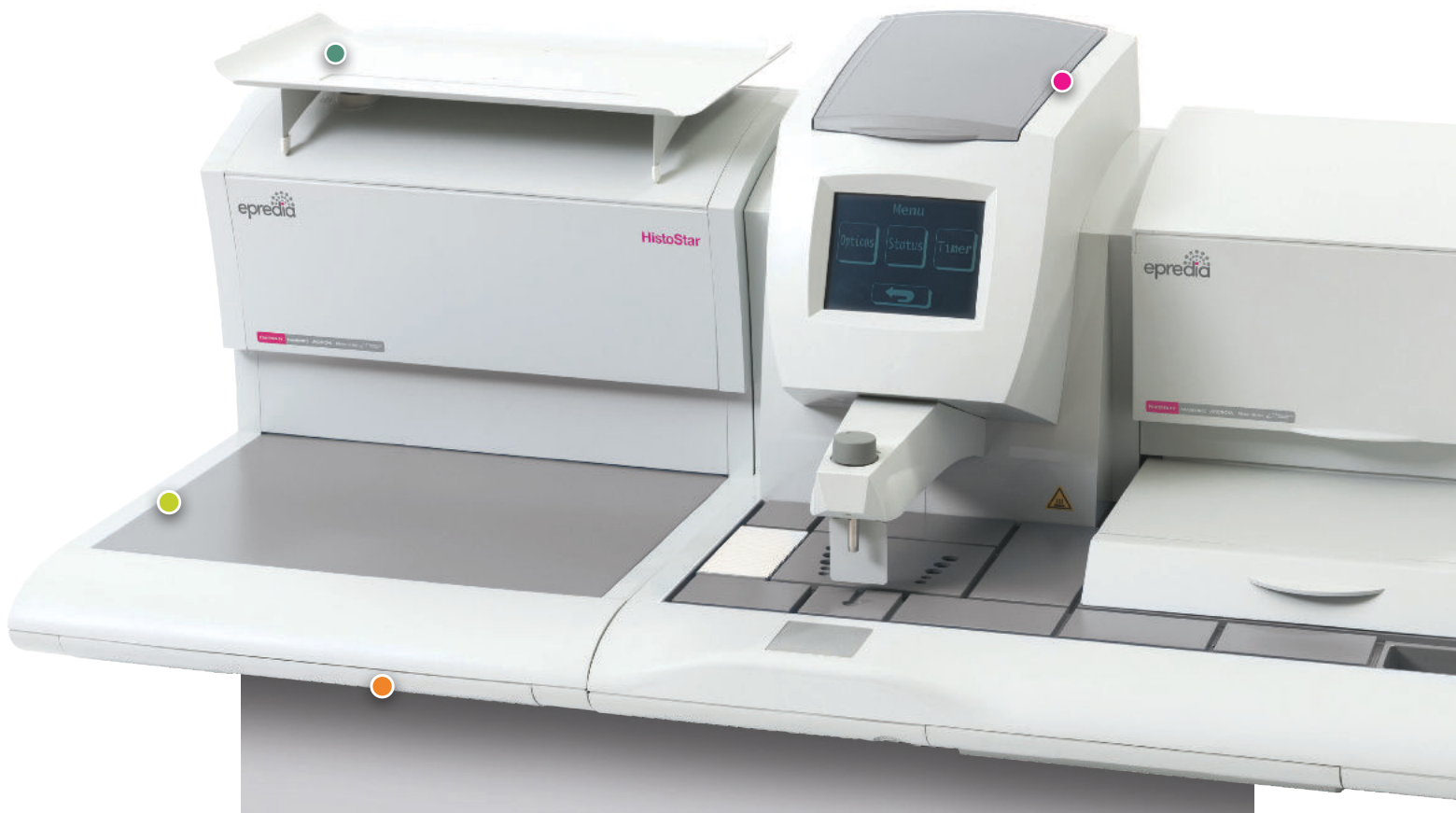

- Built-in literature tray for report storage
- Modular cold plate with 72-block capacity allows positioning to either the left or right of the embedding center for optimal workflow
- Large 5 L paraffin tank features a filtration system designed to prevent debris from entering the dispensing nozzle
- Convenient drawer for easily accessible storage

Engineered for efficiency


- Large 5 L paraffin capacity, a cold plate sized to fit 72 blocks capacity cold plate and storage for 300 cassettes.
- The Eprelia™ HistoStar™ embedding workstation is designed to minimize service requirements, ensuring dependable performance for years to come.

Flexibility to fit your workflow

- The HistoStar workstation's adjustable paraffin-dispense paddle allows users to easily adjust the paddle to the most ergonomic position.
- A variety of nozzle lengths to fit for your workflow. Even Supra Mega cassettes can be embedded with ease.
- Integrated Eprelia™ ParaTrimmer™ offers space-saving convenience. No more wax scraping at the microtome.
- The system's efficient wax residue pathway helps maintain a clean workspace.

Cool comfort

- Smooth, contoured form factor to eliminate uncomfortable pressure points.
- All user contact points are insulated to provide all-day comfort with a cool workspace.
- LED lighting uniformly illuminates the workspace reduces fatigue and minimizes errors without the clutter of awkward remote lamps. Lighting intensity is easily adjustable with five levels of illumination that can be individually controlled for the specimen and accessory areas.

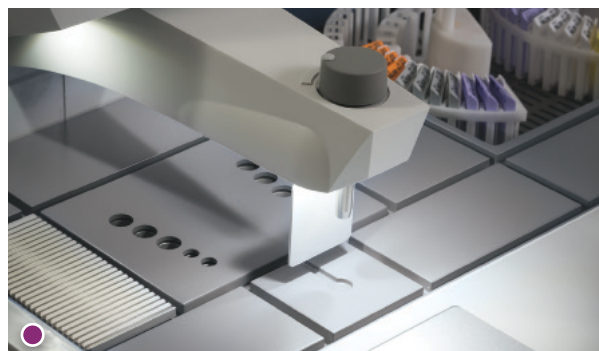

- User-adjustable cold plate temperature designed to reduce the occurrence of cracked blocks
- High-capacity storage for 200 standard base molds, 300 cassettes, 6 Eprelia Excelsior™ or 4 Revos™ baskets
- LED lighting with five brightness settings offers optimal illumination for both the specimen and accessory area
- Accessory area offers ten tool positions with five different hole sizes
- ParaTrimmer design, optimized for maximum productivity
- Intuitive, multi-lingual, touch screen control system

Epredia HistoStar Embedding Workstation Specifications

Dimensions	Embedding module (W, D, H)	25.6 x 23.6 x 16.1 in (650 x 600 x 410 mm)
	Cold plate (W, D, H)	17 x 23.6 x 16.1 in (430 x 600 x 410 mm)
Weight	Embedding module	55.12 lbs (25 kg)
	Cold plate	44.10 lbs (20 kg)
Operating noise level		<65 dB
Capacity, wax reservoir		1.32 gallon (5 litres)
Temperatures	Wax reservoir	122 - 158 °F (50 - 70 °C)
	Cold spot	41 °F (5 °C)
	Hot spot	122 - 158 °F (50 - 70 °C)
	Tissue storage	122 - 158 °F (50 - 70 °C)
	Mold storage	122 - 158 °F (50 - 70 °C)
	Cold plate	Adjustable, 10 - 27 °F (-12 to -3 °C) with a tolerance of ± 1.8 °F (1 °C)



Epredia HistoStar
Embedding Workstation

Ordering Information

Description	Order Number	
HistoStar dispensing unit with cold plate	220 - 240 Vac 50/60Hz	A81000001
	100 - 120 Vac 50/60Hz	A81000002
	220 - 240 Vac 50/60Hz (China Only)	A81000003
HistoStar cold plate	220 - 240 Vac 50/60Hz	A81000004
	100 - 120 Vac 50/60Hz	A81000005
	220 - 240Vac 50/60Hz (China Only)	A81000006
HistoStar embedding module	110 - 120 Vac / 220 - 240 Vac 50/60Hz	A81000007
	220 - 240 Vac 50/60Hz (China Only)	A81000008
Tamping tools kit, set of 2: 1 X 20 mm and 1 X 12 mm		A81010030
Spatula kit, set of two		A81010031
Cassette lid storage cover		A81010026
Magnifier assembly		A81010032
Footswitch		AP16629
Electric, heated forceps, red, 1 mm tip		A81010033
Electric, heated forceps, yellow, 2 mm tip		A81010034
Electric, heated forceps, blue, 4 mm tip		A81010035
Document holder kit		A81010036

Efficient embedding requires an outstanding paraffin

Delivery, reliability, unmatched technical support and superior product quality are all part of what you receive with every case of Epredia™ HistoStar™ paraffin. Made of paraffin and plastic polymers, these three paraffins deliver excellent processing and yield thin, flat, wrinkle-free sections.

Contact your local Epredia representative for a side-by-side comparison.



Description	Quantity	Order Number
HistoStar PE	8 x 1 kg bags	8330
HistoStar IM	8 x 1 kg bags	8331
HistoStar LP	8 x 1 kg bags	8332

Find out more at [epredia.com](https://www.epredia.com)

epredia Enhancing precision
cancer diagnostics