



# Efficient, Versatile, Ergonomic

## Epredia Histostar Embedding Center



Distributed by:  
**Axlab A/S**  
+45 35 43 18 81  
axlab@axlab.dk  
www.axlab.dk

**Axlab**

YOUR VISION.  
OUR INSPIRATION

  
**epredia**

# Epredia HistoStar Embedding Workstation

Cool, comfortable and in control



- Built-in literature tray for report storage
- Modular cold plate with 72-block capacity allows positioning to either the left or right of the embedding center for optimal workflow and features an antimicrobial coating
- Large 5 L paraffin tank features a filtration system designed to prevent debris from entering the dispensing nozzle
- Convenient drawer for easily accessible storage

## Engineered for efficiency


- Large 5 L paraffin capacity, a cold plate sized to fit 72 blocks capacity cold plate and storage for 300 cassettes.
- The Epredia™ HistoStar™ embedding workstation is designed to minimize service requirements, ensuring dependable performance for years to come.

## Flexibility to fit your workflow

- The HistoStar workstation's adjustable paraffin-dispense paddle allows users to adjust the paddle to the most ergonomic position easily.
- A variety of nozzle lengths to fit your workflow, including the Supa Mega cassette.
- Integrated Epredia™ ParaTrimmer™ offers space-saving convenience.
- Maintaining a clean workspace with an efficient wax residue pathway.

## Cool comfort

- Smooth, contoured design, eliminating uncomfortable pressure points, and insulated to provide all-day comfort with a cool workspace.
- Lighting intensity is easily adjustable with five levels of LED illumination for the specimen and accessory areas, uniformly illuminating your workspace, reducing fatigue, and minimizing errors.

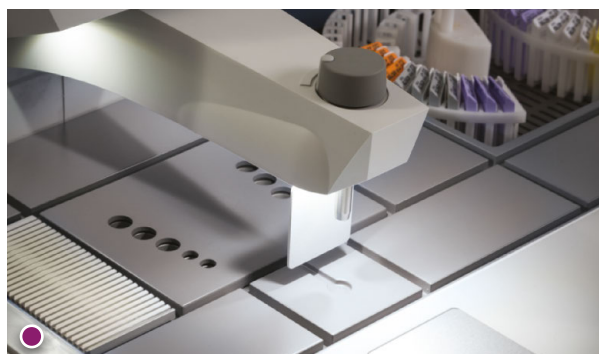

- User-adjustable cold plate temperature designed to reduce the occurrence of cracked blocks
- High-capacity storage for 300 standard base molds, 300 cassettes, and 6 Epredia Revos™ or Excelsior™ baskets
- LED lighting with five brightness settings offers optimal illumination for both the specimen and accessory area
- Accessory area offers ten tool positions with five different hole sizes
- ParaTrimmer design, optimized for maximum productivity
- Intuitive, multi-lingual, touch screen control system

## Epredia HistoStar Embedding Workstation Specifications

Dimensions	Embedding module (W, D, H)	26 x 26 x 17 in (650 x 645 x 435 mm)
	Cold plate (W, D, H)	17 x 26 x 15 in (430 x 645 x 380 mm)
Weight	Embedding module	55.12 lbs (25 kg)
	Cold plate	44.10 lbs (20 kg)
Operating noise level		>65 dB
Capacity, wax reservoir		1.32 gallon (5 litres)
Temperatures	Wax reservoir	122 - 158 °F (50 - 70 °C)
	Cold spot	41 °F (5 °C)
	Hot spot	122 - 158 °F (50 - 70 °C)
	Tissue storage	122 - 158 °F (50 - 70 °C)
	Mold storage	122 - 158 °F (50 - 70 °C)
	Cold plate	Adjustable, 10 - 27 °F (-12 to -3 °C) with a tolerance of $\pm 1.8$ °F (1 °C)



Epredia HistoStar  
Embedding Workstation

## Ordering Information

Description		Order Number
HistoStar dispensing unit with cold plate	220 - 240 Vac 50/60Hz	A81000001
	100 - 120 Vac 50/60Hz	A81000002
	220 - 240 Vac 50/60Hz (China Only)	A81000003
HistoStar cold plate	220 - 240 Vac 50/60Hz	A81000004
	100 - 120 Vac 50/60Hz	A81000005
	220 - 240Vac 50/60Hz (China Only)	A81000006
HistoStar embedding module	110 - 120 Vac / 220 - 240 Vac 50/60Hz	A81000007
	220 - 240 Vac 50/60Hz (China Only)	A81000008
Tamping tools kit, set of 2: 1 X 20 mm and 1 X 12 mm		A81010030
Spatula kit, set of two		A81010031
Cassette lid storage cover		A81010026
Magnifier assembly		A81010032
Footswitch		AP16629
Electric, heated forceps, red, 1 mm tip		A81010033
Electric, heated forceps, yellow, 2 mm tip		A81010034
Electric, heated forceps, blue, 4 mm tip		A81010035
Document holder kit		A81010036

## Efficient embedding requires an outstanding paraffin

Delivery, reliability, unmatched technical support and superior product quality are all part of what you receive with every case of Epredia™ HistoStar™ paraffin. Made of paraffin and plastic polymers, these three paraffins deliver excellent processing and yield thin, flat, wrinkle-free sections. Contact your local Epredia representative for a side-by-side comparison.

Description	Quantity	Order Number
HistoStar PE	8 x 1 kg bags	8330
HistoStar IM	8 x 1 kg bags	8331
HistoStar LP	8 x 1 kg bags	8332



Find out more at [epredia.com](https://www.epredia.com)